

Ordering Information

This product is part of the complete line of wound care products by Hollister which includes wound care dressings, wound drainage collectors, tube attachment devices and accessory products.

Educational videos and written materials are available from Hollister to help care givers learn the basics of wound care.

Wound Drainage Collectors - Sterile

Wounds up to	Code
7.6.cm	9701
9.5cm	9702
10 x 20cm	9703

Wound Drainage Collectors - Non Sterile

Wounds up to	Code
7.6.cm	9773
9.5cm	9776
10 x 20cm	9778

Wound Drainage Collectors



*Solving Draining
Wound and Fistulae
Management Problems*



For further information, please call:
FreeCall 1800 335 911 (Australia)
FreeCall 0800 167 866 (New Zealand)

Check out our Website @ www.hollister.com.au
 or www.ostomyworld.com



Australia
 6/345 Ingles Street
 Port Melbourne, Victoria 3207
New Zealand
 58 Richard Pearse Drive
 Airport Oaks, Auckland
www.ostomyworld.com.au





Wound Drainage Collectors are wound dressings which are designed to effectively collect and contain exudate from draining wounds and fistulae. The collectors are made of transparent, Odour-Barrier film. Each collector features an attached Access Cap for easy inspection and cleaning of the wound. The Convenience Drain allows fluids in the collector to be emptied and measured. The attached flexible Synthetic Skin Barrier provides long-wearing skin protection and adheres securely to the skin to prevent leakage.

Wound Drainage Collectors are Designed to Meet Clinicians' Needs

Exclusive cross-linked hydrocolloid skin barrier has demonstrated excellent clinical performance:

- Absorbs exudate.
- Conforms to hard to manage body contours; cut to fit to the exact wound shape and size.
- Adheres to skin for up to 7 days.

Patient Comfort

- Keeps the patient dry by collecting exudate; eliminates leakage between pouch and skin barrier.
- Eliminates odour by virtue of the Odour-Barrier film.
- Soft, flexible collector conforms to the patient's body.
- Skin barrier resists corrosive drainage and protects the periwound skin.

Easy to Use, Convenient and Cost Effective

- Easy to apply; virtually no residue remains on the periwound skin upon removal.
- Skin barrier can be cut to fit to the exact wound shape and size.
- Easy access is provided to inspect or treat the wound through the hinged Access Cap.
- Exudate is easily measured and/or drained through the Convenience Drain.
- Flexible skin barrier adheres to skin for up to 7 days, eliminating frequent changes.

Clinical Indications

- Sterile collectors are indicated on post-operative wounds and fistulae where exudate requires containment.
- Reduces exposure of the wound to bacteria and contaminants.
- Provides visibility and access to the wound through the Access Cap without the need to remove the collector.



Small
Accommodates wounds up to 7.6cm



Medium
Accommodates wounds up to 9.5cm



Large
Accommodates wounds up to 10cm x 20cm




Medium and large collectors are designed with gusseted sides to accommodate greater volume of exudate.



Convenience Drain can be removed with scissors to accommodate thick drainage. The collector can be closed with a Hollister Drainable Pouch Clamp. The drain will also accept tubing for continuous draining.



Universal Catheter Access Port allows the passage of a drainage catheter through the collector and provides a secure seal between the catheter and the collector.

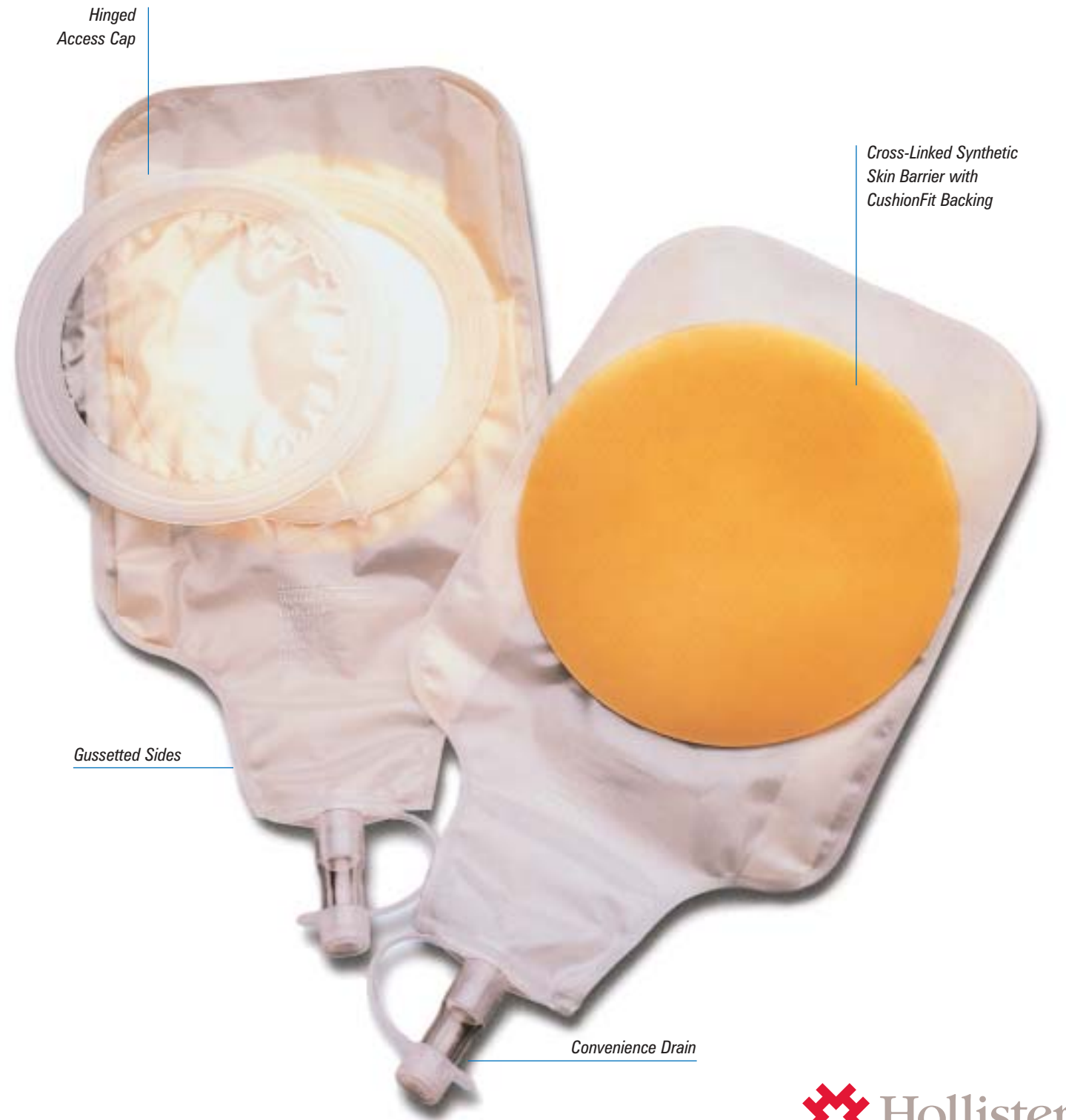

Hinged Access Cap

Gusseted Sides

Convenience Drain

Cross-Linked Synthetic Skin Barrier with CushionFit Backing