



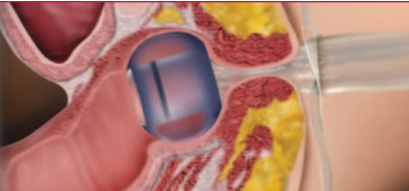

















Zassi

Bowel Management System

Unique Features

	Indication/Feature	Importance	Zassi Bowel Management System
	Indicated for Colonic Irrigation Access	<ul style="list-style-type: none"> • Provides access to the bowel for colonic irrigation 	✓
	Used to Administer Medications/Enema	<ul style="list-style-type: none"> • Permits the delivery and retention of rectally administered medications 	
	Catheter Collapse-Resistant Annulus	<ul style="list-style-type: none"> • Helps reduce catheter occlusion • Designed to help minimize leakage and expulsion 	✓
	Intraluminal Balloon	<ul style="list-style-type: none"> • Inflated for use as an introducer tip • Deflates to allow flow of stool • Inflates to occlude flow for medication and irrigant retention 	✓
	Radiopaque Marker	<ul style="list-style-type: none"> • Radiopaque Marker easily identified under X-ray or fluoroscopy 	✓
	Low-Pressure Retention Cuff	<ul style="list-style-type: none"> • Assists in holding tube in place in rectal vault 	✓ Sphere-Shaped
	Transsphincteric Zone	<ul style="list-style-type: none"> • Soft collapsible material • Designed to help avoid impact on rectal sphincter • Available in two sizes (4cm and 6cm) 	✓
	Anchor Straps	<ul style="list-style-type: none"> • Stabilizes catheter to help reduce internal migration • Helps prevent twisting or inadvertent dislodgement 	✓

	Indication/Feature	Importance	Zassi Bowel Management System
	Triple Lumen Connectors	<ul style="list-style-type: none"> • Color-coded and clearly labeled • Accepts standard luer tip syringes • Pilot balloons indicate inflation status 	✓
	Stop Flow Connector	<ul style="list-style-type: none"> • Used to inflate with air to fill the intraluminal balloon 	✓
	Irrigation/Rx Connector	<ul style="list-style-type: none"> • For administration of irrigation or medication • Adapter detaches for use with catheter tip syringe or non-IV compatible irrigation tubing 	✓
	Retention Cuff Connector	<ul style="list-style-type: none"> • Used to inflate with water to fill retention cuff 	✓
	Parylene-Coated Tube	<ul style="list-style-type: none"> • Unique coating inside and outside of tube • Helps reduce friction, and promotes drainage and odor control 	✓
	Flush/Sampling Port	<ul style="list-style-type: none"> • Provides easy access to flush drainage tubing • Split septum for convenient sample collection 	✓
	Adjustable Sheet Clip	<ul style="list-style-type: none"> • Secures tube in desired position 	✓
	Twist/Lock Connector	<ul style="list-style-type: none"> • Securely connects to either a drainable or closed collection container 	✓
	Drainable Collection Bag with Tubing and Bag Caps	<ul style="list-style-type: none"> • Can be used over a multi-day period; designed to help prevent leakage during emptying 	✓ 3000cc Capacity

	Indication/Feature	Importance	Zassi Bowel Management System
	Single-Use Collection Bag	<ul style="list-style-type: none"> • For use when a closed system is desired 	 2000cc Capacity
	Irrigation Bag	<ul style="list-style-type: none"> • For gravity-flow irrigation 	 1000cc Capacity

Note: Federal (USA) law restricts this device to sale by or on the order of a physician. Prior to use of the Zassi Bowel Management System, be sure to read: (i) the entire Zassi Bowel Management System Instructions for Use package insert supplied with the product for device Intended Use, Description, Contraindications, Warnings, Precautions, Adverse Events, and Instructions For Use; and (ii) all other package inserts and labels supplied with the product and accessories.



Attention to Detail. Attention to Life.

Hollister Incorporated
 2000 Hollister Drive
 Libertyville, Illinois 60048 USA
 1.800.323.4060

Distributed in Canada by
Hollister Limited
 95 Mary Street
 Aurora, Ontario L4G 1G3
 1.800.263.7400
www.hollister.com